
# Sefam Analyze

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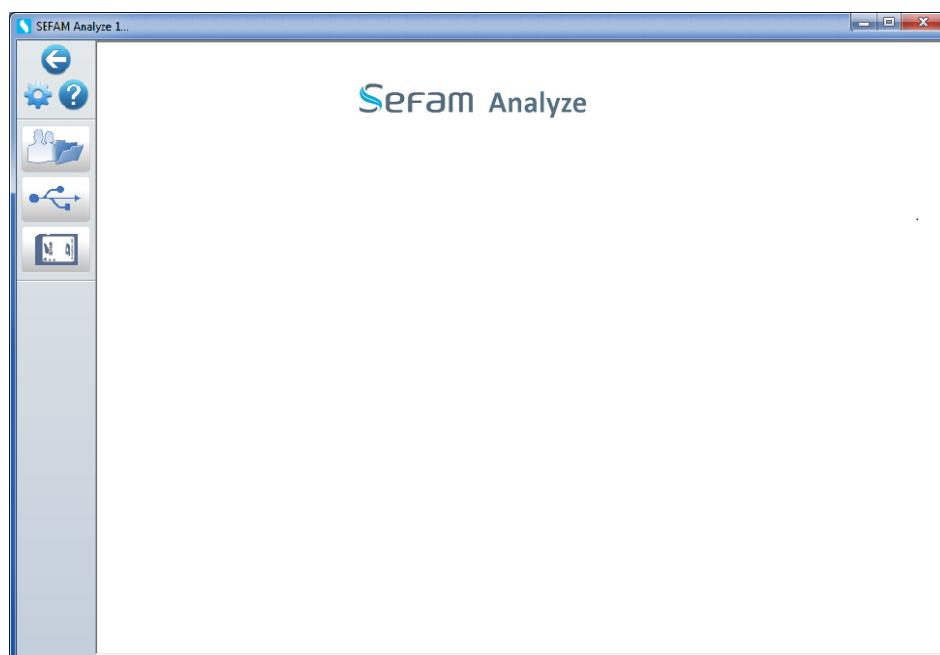# Installing the SEFAM Analyze software

## Introduction

This procedure allows to install the version 1.x.0 of the SEFAM Analyze software onto a computer equipped with a USB port or an access to the website [www.sefam-medical.com](http://www.sefam-medical.com).

## User manual

The User manual is included in the software and can be accessed by selecting the help file embedded in the software. A paper copy of the user manual in the language of your country (P/N M-156DFU01-40) may be obtained quickly and free of charge by request to your provider.

## Storing and using the USB key

Take the following precautions when storing and using the USB key.

- Check packaging integrity.
- Keep well away from humid conditions
- Clean with a clean, soft cloth.

## Installation requirements

SEFAM Analyze is 32-bit software designed for use on PCs and compatible computers.

The minimum configuration required for your PC is as follows:

- Operating system Windows XP SP3 (PII 350MHz at least 128MB)
- Windows Vista 32 bits. Home Basic ( 32 bit processor/ 512 MB of system memory )
- Windows Vista 32 bits Premium/Business/Ultimate (32 bits processor / 1 GB of system memory)
- Windows 7 32/64 bits Home Premium/ Business/Ultimate (1 GHz processor /1GB RAM, 16GB disk space)
- Windows 8, 8.1 and Windows 10
- USB port
- SD card reader.

The software can communicate via Bluetooth with the machine if the computer has this ability. The minimum supported version is Microsoft Windows 7.

## Installing the software

Before installing the software, confirm that your computer configuration meets the minimum installation requirements for the SEFAM Analyze software as defined under "Installation requirements".

### Notes:

- If the software is to be installed on a network or in an environment involving user profiles with limited access rights, it can only be installed by the network administrator who alone has the access rights necessary to correctly install the software.
- It is strongly recommended that you close all open Windows applications before proceeding with installation.
- It is recommended that you disable all screen savers, hibernation and power management on your computer.
- In case of reinstall operation, it is strongly recommended to make a copy of your patient files directory.

The SEFAM Analyze software is installed as follows:

1. Insert the USB key into the corresponding port, then open the disk and run the installation file manually.

The SEFAM Analyze installation program will start automatically. Then you can choose the language of installation in the drop down list proposed.

**Notes:**



If you install the software from the website: you must extract the file, then run the installation file manually from the directory in which the extraction was performed.

2. Follow the instructions provided in the Welcome window.

By default, the software proposes to install SEFAM Analyze in the C:\SFA directory of your hard drive (C:). If you wish to install the software in a different directory, click on Browse and select the directory. Once you have selected the directory, click on Next.

3. Then you can select on the screen the items you want to install by checking for each of them the box on the left.




- SEFAM Drivers: will install device drivers
- OMNIKEY Drivers: will install drivers for OMNIKEY card reader
- Software Options :
  - ♦ Set PDF password: options of PDF reports can be chosen
  - ♦ Size and weight units: units of weight and size can be selected
  - ♦ Analysis parameters: values of analysis parameters can be selected.
  - ♦ Setup on multi-computers: possibility to memorize previous parameters for a multi-computers deployment. The multi-computers deployment is possible if the setup directory is the same for successive installations (it is advised to make a copy on a USB key which will be used for installations).

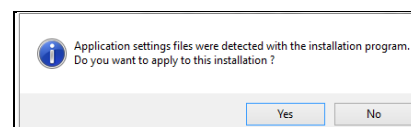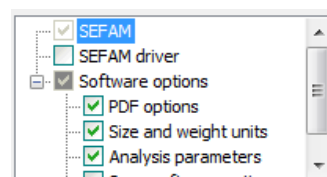**Note:**

If a multi-computers setup is detected, a dialogue box appears. Simply answer No to the question, and you will come back to the previous window.

4. Depending on what you have selected, the drivers are successively installed, then you can enable the password protection of the pdf exported document (a.) and/or select the size and weight units (b.) and/or the analysis parameters (c.) in the dialogue boxes displayed chronologically:

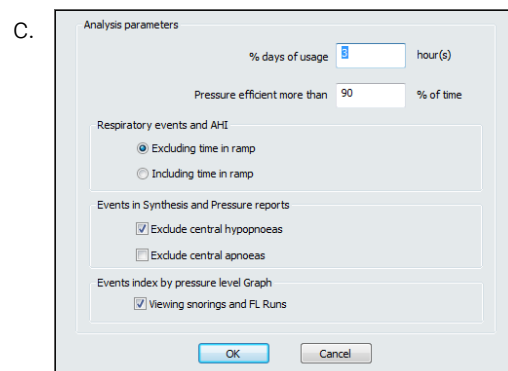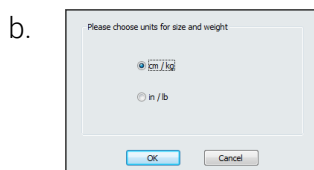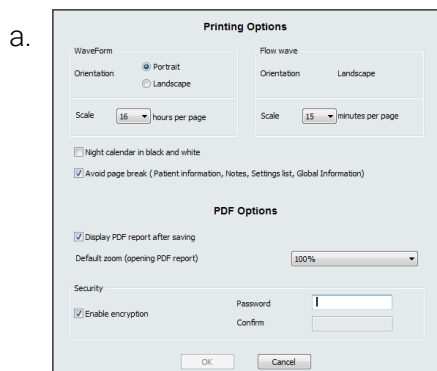5. A message informs you that the installation is completed and asks if you want to run the software. Click Yes in this case.
6. The software is installed in a folder named "SFA" by default, but you may rename it if you want.
7. When the installation is finished, a file will appear on your Desktop. It contains the program shortcut which can be accessed by double-clicking on the file.

8. This software requires Acrobat Reader to produce reports. An English version is available in the \ADOBE directory of the USB key. To install it, launch the AdbeRdr810\_en\_US.exe program file and follow the instructions for installation.

**Notes:**

- For the first installation, you must also install the drivers for the memory card reader. Follow the installation instructions described in "Driver Installation Guide" provided on the USB key, (file "Drivers installation guide.pdf").
- The drivers need to be installed only one time on the user's PC.
- Keep the USB key in the drive until the end of installation process.
- You need administrative user's rights to install the drivers.

## Installation of the PDF document viewer

This software requires Acrobat Reader to visualize reports. An English version is available in the \ADOBE directory of the USB key. To install it, launch the AdbeRdr810\_en\_US.exe program file and follow the instructions for installation.

## Uninstalling the software

**Note:**

It is strongly recommended to make a copy of your patient files directory before any uninstall operation.

To uninstall the software from the computer on which it is installed, it is recommended that you use the operating system's uninstall feature (under Control Panel, click on Add/Remove Programs and look for the SEFAM Analyze program) or go to the install directory and click on icon "uninst\_sfa.exe".

**Note:**

The uninstall procedure does not delete any SEFAM patient files on the computer's hard drive. Therefore, if you want to completely uninstall the software, it is recommended that you manually remove it from Windows Explorer by deleting the install directory.

## Upgrading the software

**Note:**

It is strongly recommended to make a copy of your patient files directory before any upgrade operation.

To upgrade the software from a previous version, you need to uninstall the previous version from the Control Panel >> Add/remove Programs, and then to run the steps 1 to 6 as described before in § "Installing the software".

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