

PolyLink

A unique integrated polygraphy system for S.Box

T R E A T M E N T






- ▷ Synchronization of signals with data from S.Box
- ▷ Automatic start
- ▷ Wireless sensors
- ▷ Sefam Analyze software

PolyLink is an accessory of S.Box. PolyLink records all the necessary parameters for the realization of a control polygraphy under treatment.

Operation

PolyLink integrates the sensors for the thoracoabdominal movement, body position and SpO₂ measurements. It ensures the acquisition and synchronization of these signals in real time with the information recorded by S.Box (flow rate, leaks...).

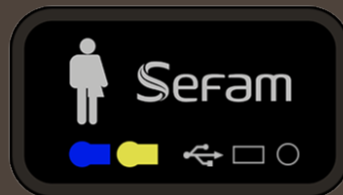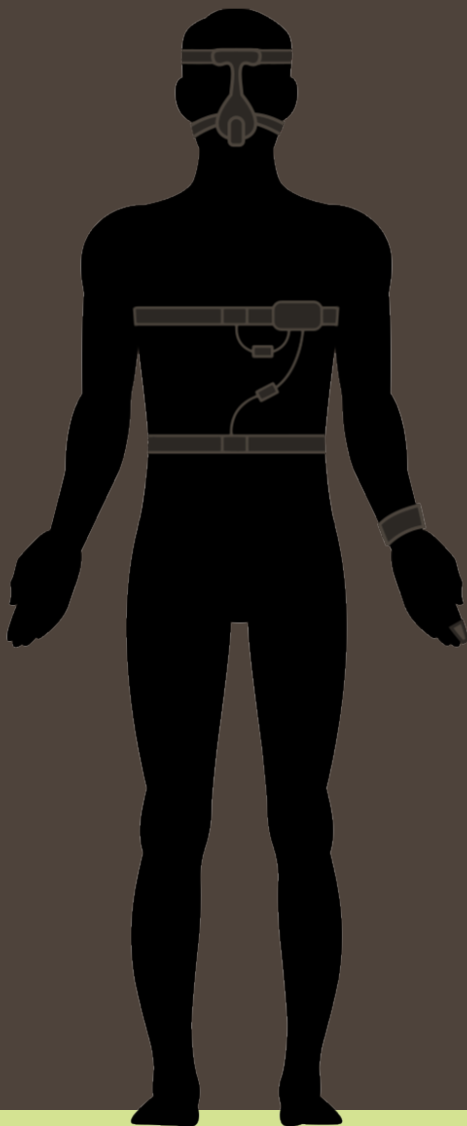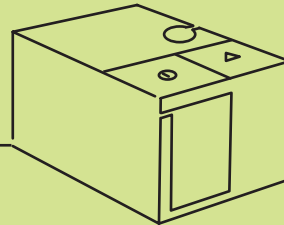

Respiratory effort


Body position





SpO₂

Heart rate

Pulse wave

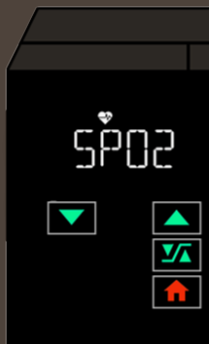


Thanks to the S.Box touchscreen users can check the proper fit of the sensors and confirm the detection of the components required for polygraphy as well as the correct communication with the CPAP.

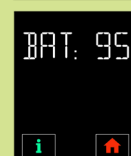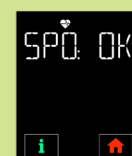

Sefam Analyze Software

As of version 2.0 of Sefam Analyze you can:

- Download all raw data recorded on the S.Box SD card, including information collected by PolyLink
- Visualize the signals and automatic and manual scoring of:
 - the respiratory events (apnea and hypopnea, obstructive and central)
 - desaturations
 - other types of events...
- Edit and print polygraphy reports as well as classic CPAP analysis (pressure, leaks etc) during the same night

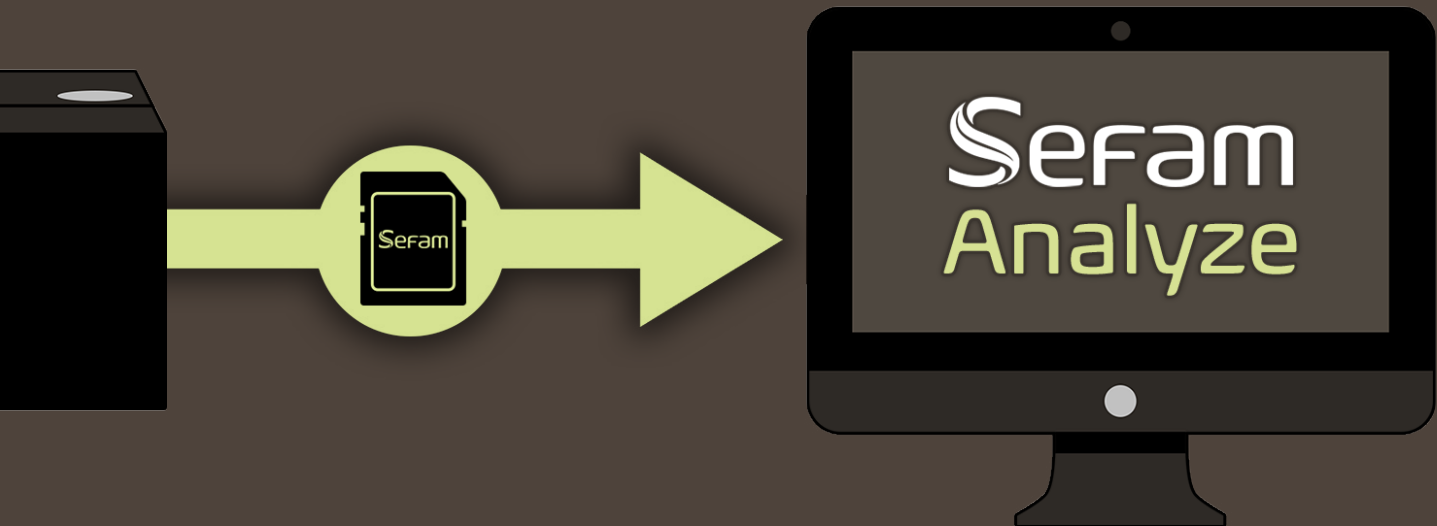

PolyLink COMPONENTS



The PolyLink effort/position kit centralizes the main data collection functions:

- ▷ connection of the THO & ABD sensors
- ▷ battery charging via USB cable
- ▷ battery power indicator
- ▷ ON/OFF switch of the system
- ▷ BLE transmission during the night

1350	POLYLINK EFFORT BELT
1368 - Kit	POLYLINK EFFORT SENSOR / ABDOMEN
1369 - Kit	POLYLINK EFFORT SENSOR / THORAX
M-116700	POLYLINK KIT (COMPLETE)
M-216730-00	POLYLINK MODULE FOR S.BOX
M-216730-01	POWER CORD FOR USB CHARGER
M-216730-02	POLYLINK EFFORT / POSITION KIT
M-216730-03	POLYLINK EFFORT / POSITION KIT HOLDER
M-416710-00	USB CHARGER FOR EFFORT / POSITION KIT
M-816705-00	POLYLINK TRANSPORT BAG



For more information:



Customer Service

T. +33 (0) 3 83 44 77 40
F. +33 (0) 3 83 44 99 87
customerservice@sefam-medical.com

www.sefam-medical.com



Head Office: SEFAM - 144 avenue Charles de Gaulle - 92200 Neuilly-sur-Seine - FRANCE
Factory and Offices: SEFAM - 10 allée Pelletier Doisy - 54600 Villers-Lès-Nancy - FRANCE



BRO208-01/05-2019

Sefam