
Sefam Analyze

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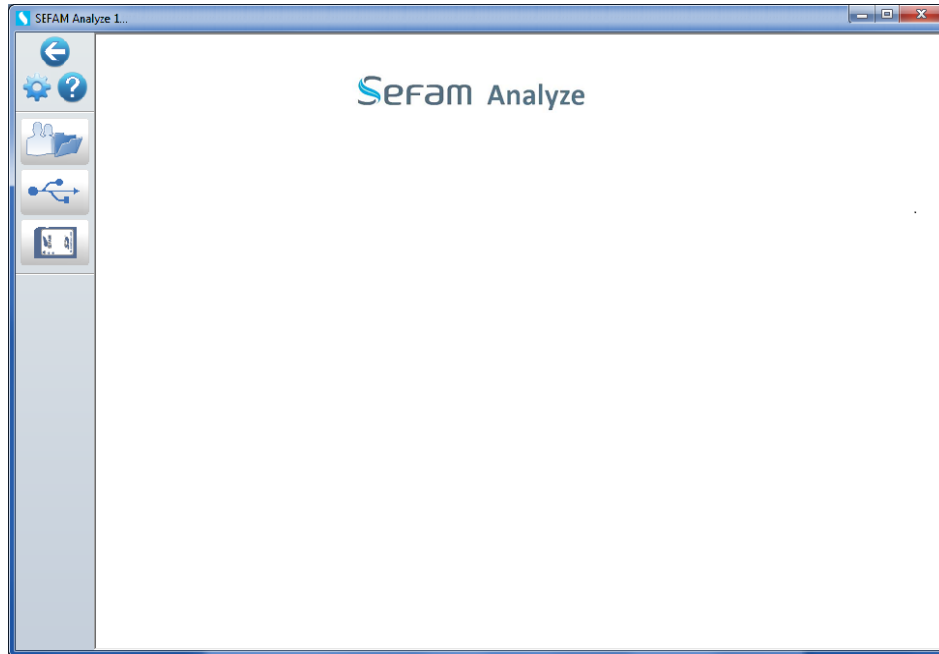Installing the SEFAM Analyze software

Introduction

This procedure allows to install the version of the SEFAM Analyze software onto a computer equipped with a USB port or an access to the website www.sefam-medical.com.

User manual

The User manual is included in the software and can be accessed by selecting the help file embedded in the software. A paper copy of the user manual in the language of your country (P/N M-156DFU01-40) may be obtained quickly and free of charge by request to your provider.

Storing and using the USB key

Take the following precautions when storing and using the USB key.

- Check packaging integrity.
- Keep well away from humid conditions
- Clean with a clean, soft cloth.

Installation requirements

SEFAM Analyze is a 32-bit software that can be used on a computer equipped with the following minimum configuration:

- Windows 10 and Windows 11(32/64)
- USB port
- SD card reader.

The software can communicate via Bluetooth with the machine if the computer has this ability. The minimum supported version is Microsoft Windows 7.

Installing the software

Before installing the software, confirm that your computer configuration meets the minimum installation requirements for the SEFAM Analyze software as defined under "Installation requirements".

Notes:

- If the software is to be installed on a network or in an environment involving user profiles with limited access rights, it can only be installed by the network administrator who alone has the access rights necessary to correctly install the software.
- It is strongly recommended that you close all open Windows applications before proceeding with installation.
- It is recommended that you disable all screen savers, hibernation, and power management on your computer.
- In case of reinstall operation, it is strongly recommended to make a copy of your patient files directory.

The SEFAM Analyze software is installed as follows:

1. Insert the USB key into the corresponding port, then open the disk and run the installation file manually.

The SEFAM Analyze installation program will start automatically. Then you can choose the language of installation in the drop-down list proposed.

Note:


If you install the software from the website: you must extract the file, then run the installation file manually from the directory in which the extraction was performed.

2. Follow the instructions provided in the Welcome window.

By default, the software proposes to install SEFAM Analyze in the C:\SFA directory of your hard drive (C:). If you wish to install the software in a different directory, click on Browse and select the directory. Once you have selected the directory, click on Next.


3. Then you can select on the screen the items you want to install by checking for each of them the box on the left.

- SEFAM Drivers: will install device drivers
- Software options:
 - ♦ Set PDF options: options of PDF reports can be chosen
 - ♦ Size and weight units: units of weight and size can be selected
 - ♦ Analysis parameters: values of analysis parameters can be selected.
 - ♦ Setup on multi-computers: possibility to memorize previous parameters for a multi-computers deployment. The multi-computers deployment is possible if the setup directory is the same for successive installations (it is advised to make a copy on a USB key which will be used for installations).



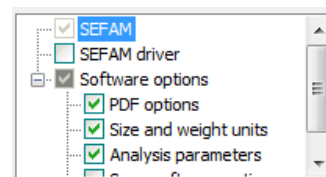
Note:

If a multi-computers setup is detected, a dialogue box appears. Simply answer No to the question, and you will come back to the previous window.

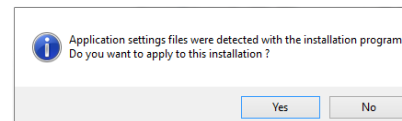4. Depending on what you have selected, the drivers are successively installed, then you can enable the password protection of the PDF exported document (a.) and/or select the size and weight units (b.) and/or the analysis parameters (c.) in the dialogue boxes displayed chronologically:

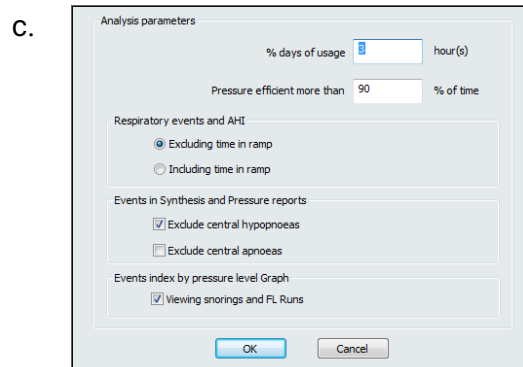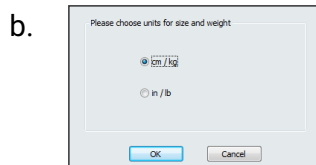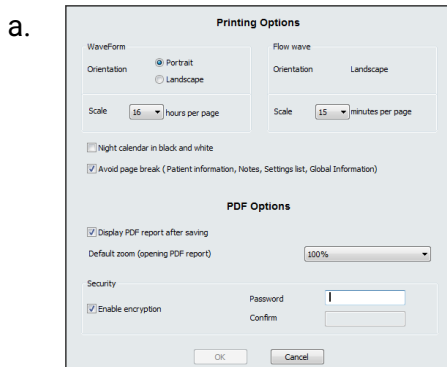5. A message informs you that the installation is completed and asks if you want to run the software. Click Yes in this case.
6. The software is installed in a folder named "SFA" by default, but you may rename it if you want.
7. When the installation is finished, a file will appear on your Desktop. It includes the program shortcut which can be accessed by double-clicking on the file.

Notes:

- For the first installation, you must also install the drivers for the memory card reader. Follow the installation instructions described in "Driver Installation Guide" provided on the USB key, (file "Drivers installation guide.pdf").
- Keep the USB key in the drive until the end of installation process.
- You need administrative user's rights to install the drivers.

Uninstalling the software

Note:

It is strongly recommended to make a copy of your patient files directory before any uninstall operation.

To uninstall the software from the computer on which it is installed, it is recommended that you use the operating system's uninstall feature (under Control Panel, click on Add/Remove Programs and look for the "SEFAM Analyze" program) or go to the install directory and click on icon "uninst_sfa.exe".

Note:

The uninstall procedure does not delete any patient files on the computer's hard drive. Therefore, if you want to completely uninstall the software, it is recommended that you manually remove it from Windows Explorer by deleting the install directory.

Upgrading the software

Note:

It is strongly recommended to make a copy of your patient files directory before any upgrade operation.

To update the software, it is advisable to uninstall the previous version (see Uninstalling the software above), then repeat steps 1 to 6 of the "Installing the software" paragraph.

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